

A Model to Study the Phenotypic Changes of Insect Cell Transfection by Copepod Super Green Fluorescent Protein (cop-GFP) in Baculovirus Expression System

Narjes Shokrollahi¹, Delavar Shahbazzadeh¹, Kamran Pooshang-Bagheri¹, Mahdi Habibi-Anbouhi², Ali Jahanian-Najafabadi³ and Mahdi Behdani^{*1}

¹Biotechnology Research Center, Venom and Biotherapeutics Molecules Lab, Pasteur Institute of Iran, Tehran, Iran;

²National Cell Bank of Iran, Pasteur Institute of Iran, Tehran, Iran; ³Department of Pharmaceutical Biotechnology, School of Pharmacy, Isfahan University of Medical Sciences and Health Services, Isfahan, Iran

Received 28 July 2015; revised 15 August 2015; accepted 26 September 2015

ABSTRACT

Background: Baculovirus expression system is one of the most attractive and powerful eukaryotic expression systems for the production of recombinant proteins. The presence of a biomarker is required to monitor transfection efficiency or protein expression levels in insect cells. **Methods:** The aim of this study was to construct a baculovirus expression vector encoding a copepod super green fluorescent protein (copGFP). In this light, the resultant vector was constructed and used for transfection of *Spodoptera frugiperda* cells. **Results:** Expression of the copGFP protein in insect cells was confirmed by fluorescent microscopy and Western-blot analysis. **Conclusion:** The application of copGFP control bacmid can be considered as an appropriate control for insect cell transfection. **DOI: 10.7508/ibj.2016.03.008**

Keywords: Baculovirus, Sf9, Transfection

Corresponding Author: Mahdi Behdani

Biotechnology Research Center, Venom and Biotherapeutics Molecules Lab, Pasteur Institute of Iran, Tehran, Iran; Tel./Fax: (+98-21) 66480780; E-mail: Behdani73042@yahoo.com

INTRODUCTION

The application of recombinant proteins has increased dramatically in the recent years^[1]. Nowadays, recombinant proteins are essential tools for studying biological processes^[2]. There is a wide variety of protein expression systems, such as bacterial, mammalian, yeast, and insect cell. Baculovirus, one of the most commonly used expression systems, is an efficient tool for high-level expression of complex eukaryotic proteins, veterinary and human vaccines as well as biopesticides^[3-5]. The *Autographa californica* nucleopolyhedrovirus is an insect-specific virus that infects predominantly the larvae of the Lepidoptera order, but usually does not infect vertebrates^[2,6]. Baculovirus expression vector is a recombinant virus genome engineered to clone the desired protein-coding sequence under the control of

strong viral promoters, polyhedron, or P10^[7]. In this study, we used Bac-to-Bac expression system that exerts Tn7 transposons for the production of recombinant baculoviruses. To monitor the transfection efficiency and expression levels of a given recombinant construct in a cell, plasmids containing reporter genes can be easily used. In the present study, we demonstrated that copepod super green fluorescent protein (copGFP), as a non-invasive, sensitive and convenient selectable marker, can be applied for direct measurement and detection of gene transfer efficiency in insect cells.

MATERIALS AND METHODS

Spodoptera frugiperda (Sf9) cells were grown in Grace's insect cell medium (Invitrogen, USA)

[Downloaded from ibj.pasteur.ac.ir on 2026-02-03]
[DOR: 20.1001.1.1028852.2016.20.3.1.4]
[DOI: 10.7508/ibj.2016.03.008]

Table 1. The specific and universal primer sets used in this study

Name	Sequence
Nb-BacF	5'-ACGGGATCCACAGGTGCAGCTGCAGGAGTCTGG-3'
Nb-BacR	5'-ACGCTCGAGTTATGAGGAGACGGTGACCTGG-3'
Cop-BacF	5'-ACGGGATCCGATGGAGAGCGACGAGAGC-3'
Cop-BacR	5'-ACGCTCGAGTTATTAGCGAGATCCGGTGG-3'
polh-Forward	5'-AAATGATAACCATCTCGC-3'
PFSBC- Reverse	5'-CCTCTACAAATGTGGTATGGC-3'
pUC/M13-Forward	5'-CCCAGTCACGACGTTGTAAAACG-3'
pUC/M13-Reverse	5'-AGCGGATAACAATTCACAGG-3'

supplemented with 10% fetal bovine serum and 1% penicillin-streptomycin (Gibco, USA) at 27°C.

VEGFR2-specific Nanobody (3VGR19) and copGFP genes were amplified from pHEN4-3VGR19 and pCDH-CMV-MCS-EF1-copGFP-T2A-Puro plasmids, respectively^[8]. PCR was used to amplify the genes with specific primers containing *Bam*HI and *Xho*I restriction sites (Table 1). The PCR reaction was performed under standard conditions for 35 cycles (Fig. 1A and 1B). Then, the PCR products were cloned into a pFastBacHTA donor vector. The resulting constructs, pFast-VGRNb and pFast-copGFP, were transformed into the *E. coli* strain TOP10F' (Invitrogen, USA) and confirmed by colony PCR using universal primers (polh-Forward and PFSBC-Reverse) and DNA sequencing.

The purified plasmids, pFast-VGRNb and pFast-copGFP, were transformed into *E. coli* DH10Bac. To confirm transposition, colony PCR was performed using M13 universal primers (Table 1). As shown in

Figure 1C and 1D, two bands, about 3230 bp (Bac-copGFP) and 2830 bp (Bac-VGRNb), were found after PCR amplification, and the recombinant bacmids were purified by the Carmo's method^[9].

The recombinant bacmid was transfected into insect cells using the Bac-to-Bac baculovirus expression system (Invitrogen, USA) according to the manufacturer's instructions. Briefly, 1 µg recombinant bacmid was dissolved in 100 µl Grace's insect cell culture. The resulting mixture was then mixed with 6 µl Cellfectin (Invitrogen, USA) reagent dissolved in 100 µl Grace's insect cell culture (Invitrogen, USA), after 20 min of incubation, 800 µl Grace's insect cell culture was added to the mixture, and subsequently transfected into 8×10⁵ *Sf9* cells. The transfection mixture was then removed and replaced with a fresh growth medium containing 10% FBS and 1% penicillin-streptomycin in each well. Seventy-two hours post transfection, the signs of infection were

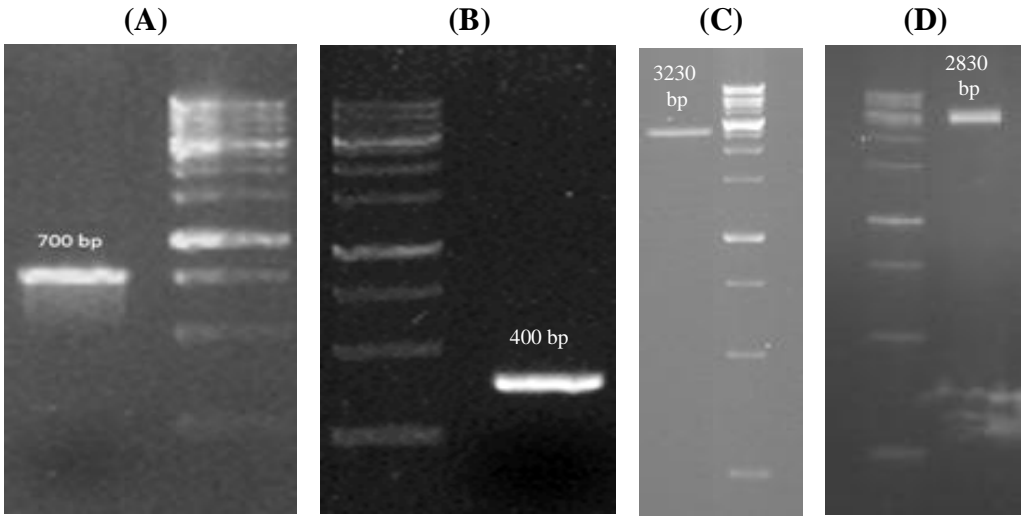


Fig. 1. PCR products analyzed on 1.5% agarose gel. As shown in the Figure, the size of copGFP (A) and 3VGR19 (B) PCR products are 700 and 400 bp, respectively. Colony PCR was performed using M13 universal primers on Bac-VGRNb (C) and Bac-copGFP (D) to confirm transposition into bacmid. The size of the DNA marker lanes from bottom to top is 250, 500, 750, 1000, 1500, 2000, 2500, 3000, 4000, 5000, 6000, 8000, and 10000 bp.

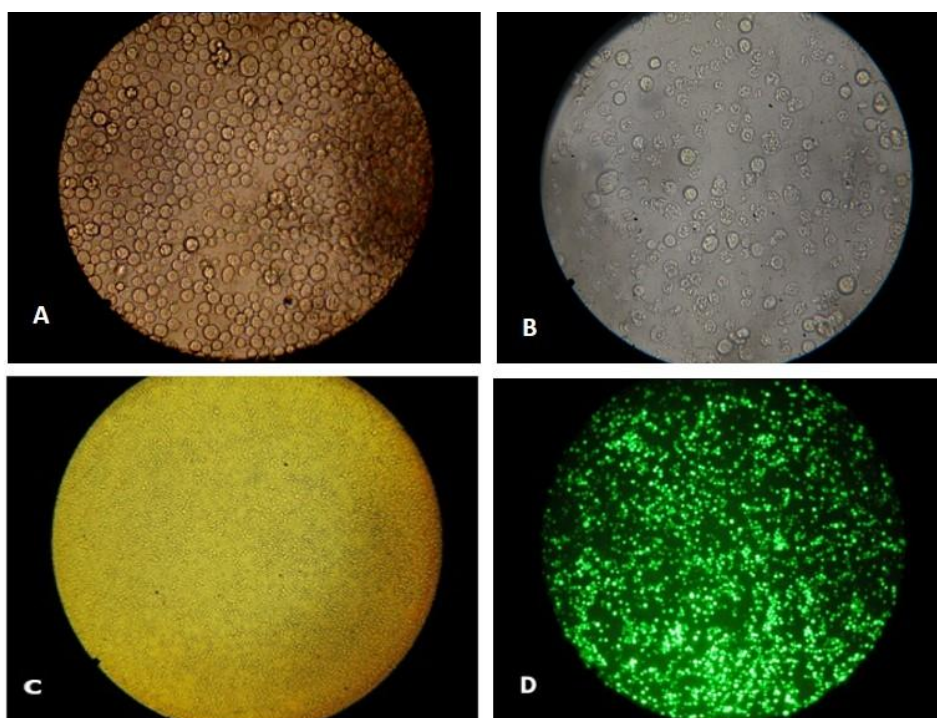


Fig. 2. Result of transfection in *Sf9* cells. A) Control cells ($\times 40$); B) Cells transfected with Bac-VGRNb ($\times 40$); C) Control cells ($\times 10$); D) Cells transfected with Bac-copGFP (fluorescent microscope) ($\times 10$).

observed under light and fluorescent microscopes. Transfected cells were analyzed after 48 h for visible hallmarks of *Autographa californica nucleopolyhedrovirus* infection. These hallmarks include cell growth arrest, cell size enlargement, and cell monolayer detachment.

A large number of transfected cells carrying Bac-copGFP were found to have a bright green color, representing the successful transfection. Cell size enlargement and cell monolayer detachment were observed in the *Sf9* cells transfected with bacmid after three days. In addition, approximately 70% of the Bac-copGFP transfected cells exhibited a bright green color (Fig. 2).

After 72 hours post transfection, the supernatant containing virus particles was recovered by centrifugation at $664 \times g$ at $4^{\circ}C$ for 3 min. The cell pellets were lysed by freeze-thaw cycles using liquid nitrogen and used for Western-blot analysis with the anti-His antibody. As shown in Figure 3, copGFP and 3VGR19 yielded 40- and 20-kDa bands, respectively.

RESULTS AND DISCUSSION

In recent years, a variety of cell expression systems has gained increasing attention for their application in recombinant protein production. However, there seems

to be a need for biomarkers to confirm the transfection process. In the case of transient transfection, there are a number of reporter genes, such as secreted alkaline phosphatase, β -galactosidase, firefly luciferase, and chloramphenicol acetyltransferase, which are used to detect and quantify transfection efficiency *in vitro*^[10]. In the above-mentioned reporter genes, the gene transfer efficiency is indirectly assessed using the conditioned medium or cellular extracts of transfected cells. In addition, *in vivo* reporter assays, such as *in situ* β -gal staining, β -glucuronidase, and luciferase, are commonly used for transfection detection in either fixed cells or tissue sections. These methods allow the direct visualization of transfected cells after staining with enzymatic substrates or antibodies. In these assays, not only the presence of fixed cell cultures but also their reaction with substrates are required to detect transfected cells^[10].

GFP, as a reporter protein, is widely used to monitor gene expression and protein localization in a broad range of cells and organisms^[10,11]. Several GFP-like proteins have been isolated from a number of copepod aquatic crustacean species^[11], the brightest of which has been made commercially available (Evrogen) under the name copGFP derived from copepod *Pontellina plumata*^[12,13]. copGFP is characterized by superbright green fluorescence (at least 1.3 times brighter than EGFP), a fast maturation rate at a wide

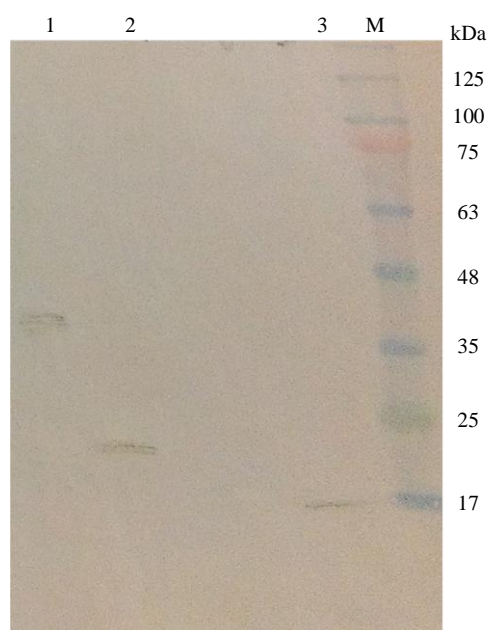


Fig. 3. Western-blot analysis for cell pellets obtained from *Sf9* cells transfected with recombinant baculovirus. Lane 1, copGFP (40 kDa); Lane 2, 3VGR19 (20 kDa); Lane 3, control protein (17kDa); M, marker

range of temperatures^[14], non-toxicity, high stability at a wide range of pH (pH 4-12), and no need for any additional cofactors or substrates^[11]. During the years, GFP has been highly regarded as a biomarker in a wide variety of cells^[15, 16].

In the current study, it was demonstrated that a reporter system containing the cop-GFP has a unique profile to monitor the transfection process in insect cells (*Sf9* cells). For this purpose, two recombinant baculovirus systems, containing copGFP as a biomarker and 3VGR19 as a control, were constructed to confirm the transfection procedure in insect cells. The present study clearly shows that the transfected cells can be easily detected using copGFP.

Following the transfection of insect cells with bacmids containing copGFP and 3VGR19, cytopathic effects were found in insect cells, and virions containing the gene of interest were released. In conclusion, the copGFP control bacmid (Bac-copGFP) can be used as a control for monitoring and optimizing transfection efficiency in insect cell transfection systems.

ACKNOWLEDGMENTS

This work was financially supported by Iran National Science Foundation (Grant No. 91004026) and Pasteur Institute of Iran (Tehran).

CONFLICT OF INTEREST. None declared.

REFERENCES

1. Terpe K. Overview of bacterial expression systems for heterologous protein production: from molecular and biochemical fundamentals to commercial systems. *Applied microbiology and biotechnology* 2006; **72**(2): 211-222.
2. Brondy WH. Selecting an appropriate method for expressing a recombinant protein. *Methods in enzymology* 2009; **463**: 131-147.
3. Nguyen Q, Chan LC, Nielsen LK, Reid S. Genome scale analysis of differential mRNA expression of *Helicoverpa zea* insect cells infected with a *H. armigera* baculovirus. *Virology* 2013; **444**(1-2): 158-170.
4. Agathos SN. Production scale insect cell culture. *Biotechnology advances* 1991; **9**(1): 51-68.
5. Feil R, Miiller S, Hofmann F. High-level expression of functional cGMP-dependent protein kinase using the baculovirus system. *FEBS letters* 1993; **336**(1): 163-167.
6. Brown T.A. Gene cloning and DNA analysis: an introduction. Wiley-Blackwell; 2010.
7. Zhu FX, Qi YP, Huang YX, Bing Q. A baculovirus vector derived from immediately early gene promoter of Autographa callifornica nuclear polyhedrosis virus. *Journal of virological methods*. 1996; **62**(1): 71-79.
8. Behdani M, Zeinali S, Khanahmad H, Karimipour M, Asadzadeh N, Azadmanesh K, Khabiri A, Schoonooghe S, Habibi Anbouhi M, Hassanzadeh-Ghassabeh G, Muyldermans S. Generation and characterization of a functional Nanobody against the vascular endothelial growth factor receptor-2; angiogenesis cell receptor. *Molecular immunology* 2012; **50**(1-2): 35-41.
9. Carmo AC, Giovanni DN, Correa TP, Martins LM, Stocco RC, Suazo CA, Moraes RH, Veiga AB, Mendonca RZ. Expression of an antiviral protein from *Lonomia obliqua* hemolymph in baculovirus/insect cell system. *Antiviral research* 2012; **94**(2): 126-130.
10. Zhang G, Gurtu V, Kain SR. An enhanced green fluorescent protein allows sensitive detection of gene transfer in mammalian cells. *Biochemical and biophysical research communications* 1996; **227**(3): 707-711.
11. Olenych SG, Claxton NS, Ottenberg GK, Davidson MW. The fluorescent protein color palette. *Current protocols in cell biology* 2007; **21**: 1-34.
12. Tennstaedt A, Aswendt M, Adamczak J, Collienne U, Selt M, Schneider G, Henn N, Schaefer C, Lagouge M, Wiedermann D, Kloppenburg P, Hoehn M. Human neural stem cell intracerebral grafts show spontaneous early neuronal differentiation after several weeks. *Biomaterials* 2015; **44**: 143-154.
13. Shagin DA, Barsova EV, Yanushevich YG, Fradkov AF, Lukyanov KA, Labas YA, Semenova TN, Ugalde JA, Meyers A, Nunez JM, Widder EA, Lukyanov SA, Matz MV. GFP-like proteins as ubiquitous metazoan superfamily: Evolution of functional features and structural complexity. *Molecular biology and evolution* 2004; **21**(5): 841-850.
14. Ro S, Park C, Jin J, Zheng H, Blair PJ, Redelman D,

- Ward SM, Yan W, Sanders KM. A model to study the phenotypic changes of interstitial cells of Cajal in gastrointestinal diseases. *Gastroenterology* 2010; **138**(3): 1068-1078.
15. Andersen JB, Sternberg C, Poulsen Lk, Bjørn SP, Givskov M, Molin S. New Unstable Variants of Green Fluorescent Protein for Studies of Transient Gene Expression in Bacteria. *Applied and environmental microbiology* 1998; **64**(6): 2240-2246.
16. Yeh E, Gustafson K, Boulianne GL. Green fluorescent protein as vital marker and reporter of gene expression in Drosophila. *Genetics* 1995; **92**(15): 7036-7040.